


AVANOS


CORFLO*
PEG Kits

WE HAVE THE RIGHT SOLUTION

For long term enteral feeding⁽¹⁾

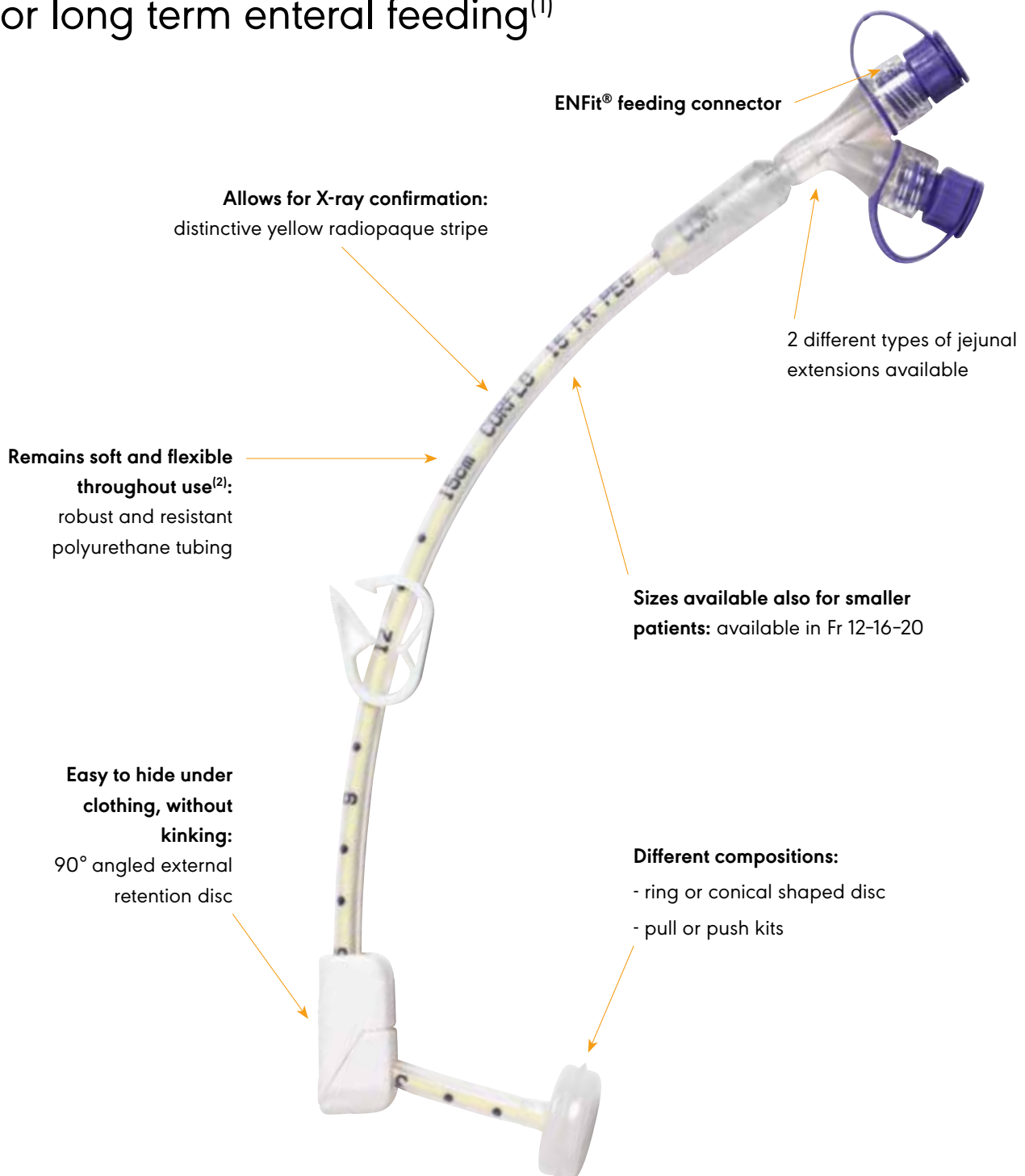

CORFLO* PEG Kits

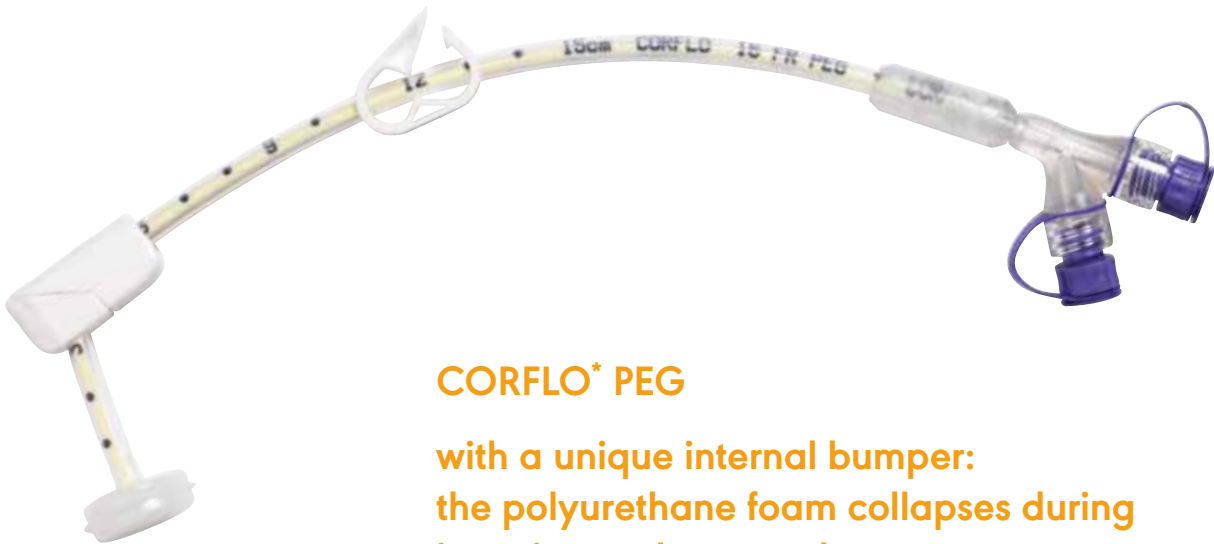

CORFLO* PEG

with a unique internal bumper:
the polyurethane foam collapses during
insertion and removal



REDUCES THE BUMPER DIAMETER SIGNIFICANTLY DURING PASSAGE THROUGH THE OESOPHAGUS:

- ideal in case of strictures e.g. esophageal atresia, cancer,...
- useful for smaller patients

ALLOWS REMOVAL BY TRACTION, WITHOUT ENDOSCOPE:

- prevents a second endoscopy and reduces potential complications arising from the anaesthesia
- minimises demands on resources

CORFLO* PEG

with a unique internal bumper:
ring shaped design with upstanding sides

REDUCES THE RISK FOR BURIED BUMPER SYNDROME⁽³⁾

especially useful when long term enteral feeding is indicated



CORFLO* PEG Additional Features:

- Polyurethane design: provides a larger inner lumen than silicone, more robust and resistant to deterioration than silicone
- Distinctive **YELLOW** radiopaque stripe tube design: allows for X-ray confirmation
- 90° angled external fixation disc holds the tube secure: prevents kinking of the tube, easy to hide under clothing
- 3 different French sizes: Fr 12-16-20
- Ring or conical internal bumper for Fr 20
- Pull or Push version



RING



CONICAL

CORFLO* PEJ Features:

- Multiple placement options for Fr 16 and Fr 20:
Grip tip design for endoscopic pull along placement.
Unique "**anti-clog tip**" with open end design for over-the-wire placement and allowing easier aspiration⁽⁴⁾.
- Centimeter markings provide a reference for placement confirmation
- Adjustable sleeve accommodates specific tube length needed for each individual without potential risk of accidental dislodgment.
- Double lumen connector:
simultaneous gastric drainage and intestinal feeding possible.
- Polyurethane provides strength and a large inner diameter.



GRIP TIP



OPEN TIP

CORFLO* PEG Kit

Each kit contains the following components:

- Polyurethane PEG tubing with traction removable internal bolster.
- Ring shaped internal retention bolster (see table for conical shaped version)
- 2 types of external retention discs
- Fast release clamp
- ENFit® feeding connector
- Towel and drape
- Scalpel
- Jelly
- Stainless steel puncture needle and trocar G18

The "E1-type" of kits have:

- An additional G14 puncture needle with teflon trocar
- No jelly

The "E2-type" of kits have:

- A blue coated pull wire
- 1 G18 puncture needle, except Fr 12 kit has a Gauge 14 puncture needle with teflon coated trocar
- No scalpel

The product item# 5020E2 contains in addition to E2-type of kits:

- 3 needles: Gauge 14, 18, 25
- Scalpel, syringe and gauzes

CORFLO* PEG kits: Pull technique – 2 eaches/packaging

SIZE	DESCRIPTION	DESIGN	CODE
Fr 12	PEG Kit Fr 12	Ring shaped disc, sterile	50-6012E1
Fr 16	PEG Kit Fr 16	Ring shaped disc, sterile	50-6016E1
Fr 20	PEG Kit Fr 20	Ring shaped disc, sterile	50-6020E1
Fr 20	PEG Kit Fr 20	Conical shaped disc, sterile	50-5020E2

CORFLO* PEG kits: Push technique – 2 eaches/packaging

SIZE	DESCRIPTION	DESIGN	CODE
Fr 16	PEG Kit Fr 16	Ring shaped disc, sterile	50-4016E1
Fr 20	PEG Kit Fr 20	Ring shaped disc, sterile	50-4020E1

CORFLO* PEG Jejunal tubes: Pull along technique – 5 eaches/packaging

SIZE	DESCRIPTION	DESIGN	CODE
Fr 6	Jejunal tube for use with CORFLO* PEG Fr 16 : Fr 6 / 91 cm, 1 gr weighted	Grip tip, non-sterile	50-7366
Fr 8	Jejunal tube for use with CORFLO* PEG Fr 20 : Fr 8 / 91 cm, 3 gr weighted	Grip tip, non-sterile	50-7368
Fr 10	Jejunal tube for use with CORFLO* PEG Fr 20 : Fr 10 / 91 cm, 3 gr weighted	Grip tip, non-sterile	50-7361

CORFLO* PEG Jejunal tubes: Over-the-wire technique – 2 eaches/packaging

SIZE	DESCRIPTION	DESIGN	CODE
Fr 8	Jejunal tube for use with CORFLO* PEG Fr 20 : Fr 8 / 109 cm, unweighted	Open tip, non-sterile	50-1438
Fr 10	Jejunal tube for use with CORFLO* PEG Fr 20 : Fr 10 / 109 cm, unweighted	Open tip, non-sterile	50-1431

CORFLO* PEG Adapter Repair Kit

SIZE	DESCRIPTION	DESIGN	CODE
Fr 12	CORFLO* PEG Adapter Repair Kit 12Fr		50-6112
Fr 16	CORFLO* PEG Adapter Repair Kit 16Fr		50-6116
Fr 20	CORFLO* PEG Adapter Repair Kit 20Fr		50-6120

References: 1. Claims based on EtQ COPY-01709. 2. Best C. Nursing Standard 2005; 20 (3): 59-65. 3. Dowman JK, et al.: Observed high incidence of buried bumper syndrome associated with Freka PEG tubes. Frontline Gastroenterology 2015; 6: 194-198. 4. Silk DBA, et al.: JPEN, 1987; 11: 378-383.

For Healthcare professional use only. Please read carefully the instructions for use provided with the Medical Device for indications, contraindications, warnings, precautions, and potential adverse effects.

AVANOS

For more information, please send an email to customerservice.uk.ie@avanos.com or visit www.avanos.co.uk.

*Registered Trademark or Trademark of Avanos Medical, Inc., or its affiliates. The COLOR YELLOW is a Registered Trademark or Trademark of Avanos Medical, Inc. or its affiliates. ©2019 AVANS. All rights reserved. ENFit® is a Registered Trademark or Trademark of Global Enteral Device Supplier Association, Inc. Used with permission. HC994-00-UK - COPY-03844